




NEW

Vertical Oval Soft Convex Pouches


The **Vertical Oval Soft Convex Pouches** have a **vertically** positioned **oval soft convex barrier**, offering more options for pouch placement than traditional convex skin barriers. Ideal for elliptical or irregular shaped **vertically** (head to toe) positioned stomas, including temporary or permanent loop ileostomies.

The benefits to you...

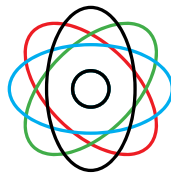

Vertical Oval Soft Convexity

- Ideal for **retracted stomas, stomas in folds and creases** and **other difficult to manage scenarios**, which don't respond to conventional **pre-cut** or **cut to fit** soft convexity.
- The **oval skin barrier** provides a more accurate fit – improved, even distribution of pressure around stoma, especially in difficult areas where skin is at risk and where leakage is more likely to occur.
- Provides gentle pressure where needed for optimum results – helping to improve skin condition, wear time, with enhanced comfort and conformity.
- Greater **flexibility** and **comfort** are provided by a 0.9mm thick skin-friendly traditional **hydrocolloid skin barrier adhesive**, backed with a **SuperFlex stretchable** foam layer providing optimum stomal protection.
- **Optimum** depth and **oval profile** for managing many irregular shaped stomas.
- Compared to traditional rigid convex skin barriers, these pouches incorporate a **softer** and more **flexible oval** insert for noticeable **comfort** when moving and bending.
- For enhanced security, each pouch incorporates a **unique 4-point belt attachment** feature for the optional use of a 50mm or 25mm wide adjustable ostomy belt.

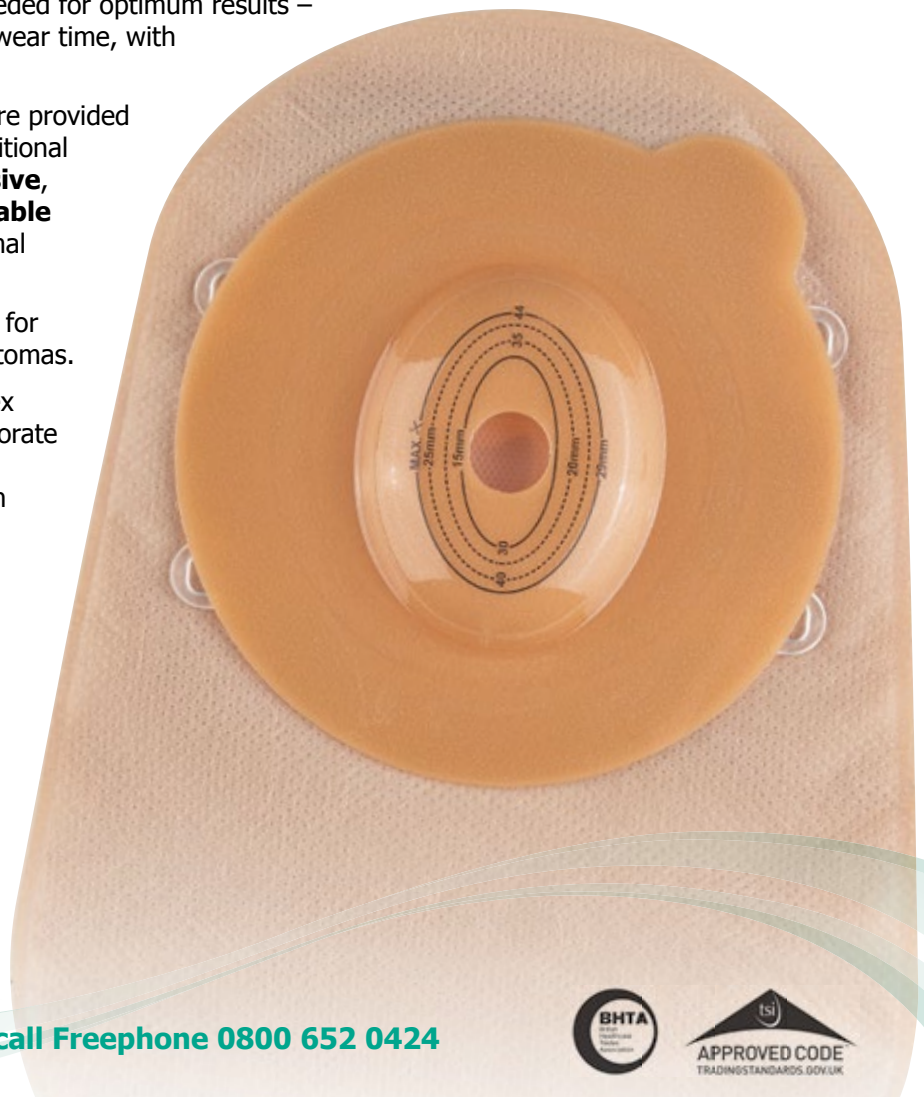

For more information please call Freephone 0800 652 0424



Ileomate Vertical Oval Soft Convex Drainable Pouches

- **Vertical oval soft convex** protective skin barrier.
- Easy **SoftSafe** closure system, providing comfort with maximum security.
- Unique **Finger Pocket** for easy drainage and hygienic cleaning of the pouch.
- Split cover to assist with pouch application and observation of the stoma.
- High capacity 32mm bar filter helps to prevent ballooning and pancaking.
- Filter protection enhanced by a **Gore-Tex®** membrane to prevent 'strike through' from liquid stool.
- Available in a large **cut to fit** pouch: **13-44 x 29mm** and an extra large **cut to fit** pouch: **13-54 x 37mm**.

Adjustable Ostomy Belts

- Available in a standard 25mm to fit all sizes or a wider 50mm belt in various waist sizes.
- Easy to attach.
- The 50mm wide belt, when attached to the **oval soft convex pouches**, will ensure even distribution of pressure and support around the peristomal area.
- The 25mm wide belt comes with a 2 to 1 adaptor to allow the use of the 4-point belt attachment on our **oval soft convex pouches**.
- Machine washable.

All oval convex products should be used upon the advice and guidance of a Stoma Care Nurse.

This product range is latex-free.

Available on NHS prescription.

Ileomate Vertical Oval Soft Convex Drainable Pouches

Ileomate

Product Codes

Large Drainable Pouch – Pack of 10

Stoma Size **Transparent**
Cut to fit 13-44 x 29mm ILTCO113-44V

Extra Large Drainable Pouch – Pack of 10

Stoma Size **Transparent**
Cut to fit 13-54 x 37mm ILTCO113-54V



Adjustable Ostomy Belts

Product Codes - Pack of 1

| Belt Width | Waist Size | Product Code |
|------------|-------------------------------------|--------------|
| 25mm | One size fits all | SB25 |
| 50mm | Small (60-75cm) | SB50S |
| 50mm | Medium (75-90cm) | SB50M |
| 50mm | Large (90-110cm) | SB50L |
| 50mm | Extra Large (110-130cm) | SB50XL |
| 50mm | Extra Extra Large (130cm and above) | SB50XXL |

